

Ver 2.
Plasma Therapy Series

PLASMATIQUE PREMIUM PLAMERE LEAF



Discover Wonders of Plasma
Illuminating Your Beauty Journey!

Plasmatique™
COLD PLASMA

PREMIUM
PLAMERE™
Esthetic Multi-Plasma Device

Leaf
FUSION PLASMA

GCS Co., Ltd.

Seoul office 8F, Magic Kingdom Bldg., 18, Teheran-ro 83-gil, Gangnam-gu, Seoul, 06169 Korea

Headquarter 3F, 33, Sagimakgol-ro 62beon-gil, Jungwon-gu, Seongnam-si, Gyeonggi-do, 13211, Korea
(Sangdaewon-dong, Center M)

Email salse@gcshcb.com

Ozone FREE!

Unmatched cold plasma experience with ozone safety technology.

Plasmatique™ makes treatments safer than other plasma devices.

It collects ozone generated during use with a unique suction mechanism to control emission and comply with safety regulations.

Choose Plasmatique™ for a safe and environmentally responsible beauty therapy.

Product Name	PLASMATIQUE™
Model Name	GPM - 1000
Output Power	MAX 7W
Input Voltage	DC 24 4A
Cold Plasma Energy Level	1.4W, 2.8W, 4.2W, 5.6W, 7W
Cold Plasma Temperature	≤55°C
Ozone Emission	Less than 0.05ppm
Size	300*260*135(Excluding protrusions)
Net weight	3Kg

GCS Co., Ltd.

Seoul office 8F, Magic Kingdom Bldg., 18, Teheran-ro 83-gil, Gangnam-gu, Seoul, 06169 Korea

Headquarter 3F, 33, Sagimakgol-ro 62beon-gil, Jungwon-gu, Seongnam-si, Gyeonggi-do, 13211, Korea
(Sangdaewon-dong, Center M)

Email sales@gcshcb.com


Ver 2.

Plasma Therapy Series

PLASMATIQUE™

Cold Plasma Device

For Skin Regeneration & Anti-Aging Therapy

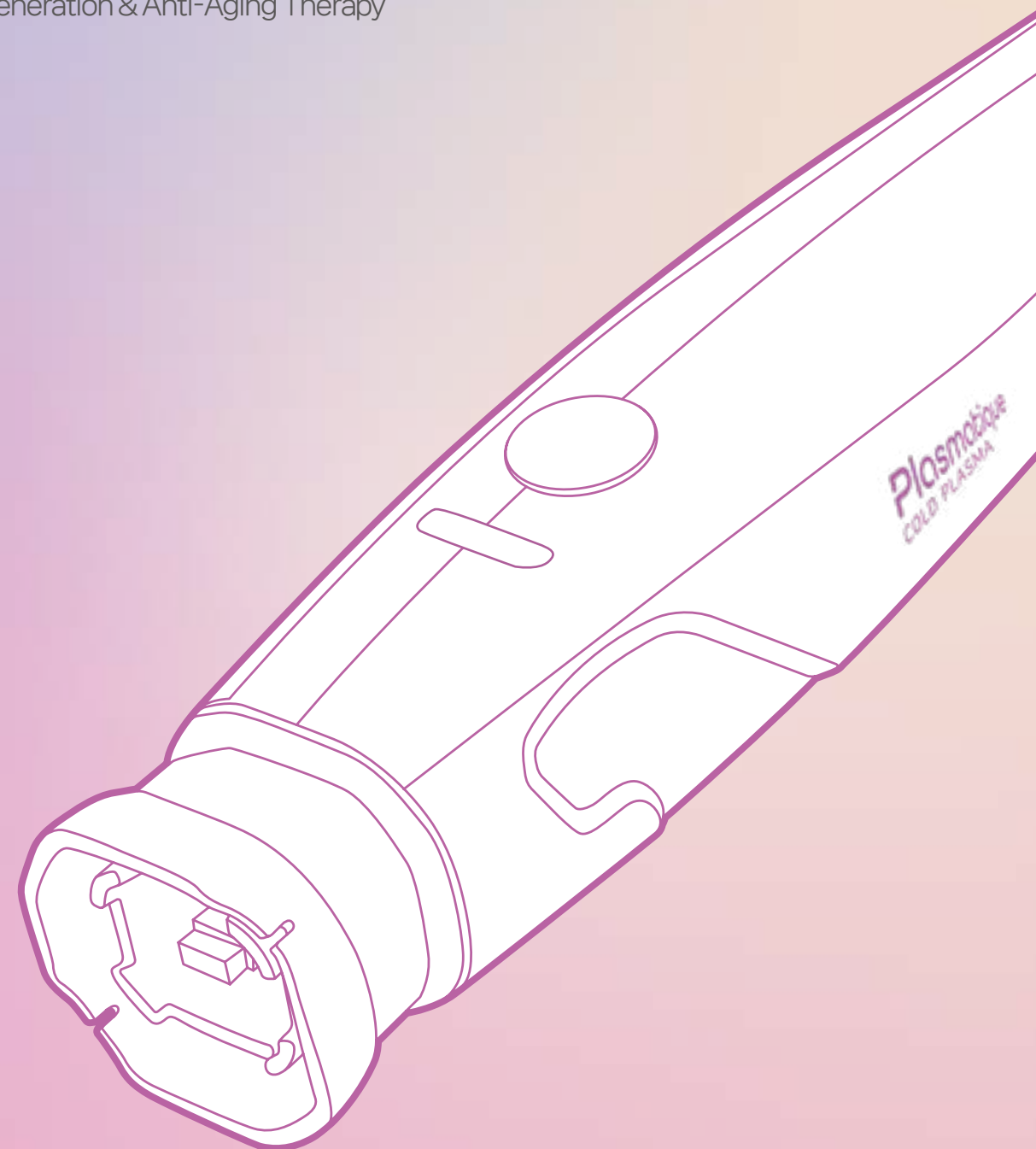

Cold Atmospheric Plasma Device for Anti-Aging

Plasmatique™ penetrates 13,000volts of Plasma into the dermal and subcutaneous layers of the skin to stimulate collagen production and remodel the cellular structure from the inside out.



01 Effects of Plasmatique™ Treatment

Plasmatique™ Boosting : Plasmatique™ temporarily disrupts cell adhesion molecules, boosting the absorption of skin boosters by up to 120 times!

Plasmatique™ Toning : Plasmatique™ works by suppressing melanin pigment, promoting a luminous, flawless complexion that radiates confidence.

Plasmatique™ Clear : Experience the soothing power of Plasmatique™'s sterilizing effect, targeting acne and atopy.


Plasmatique™ Repair : Plasmatique™ stimulates collagen and elastin synthesis, it fortifies your skin's regenerative capacity, preventing dermis layer density loss for a truly youthful glow.

02 Semi-Permanent Cold Plasma Tip

Plasmatique™ Cold Plasma tip that can be used up to 1500 hours of treatment.

03 Safe Ozone Plasma with Ozone Filter (≤0.05ppm)


FDA has set a limit on ozone by less than 0.05 ppm at any given time. Plasmatique™ is the only cold plasma device which has ozone filter and suction in the world. You can enjoy a safe cold plasma treatment without "Ozone" smell.

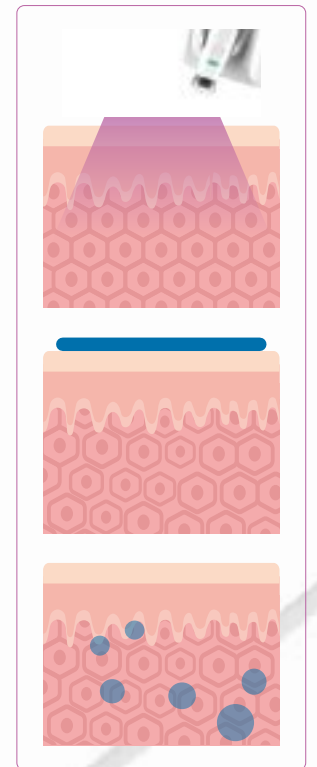

Ver 2.
Plasma Therapy Series

PREMIUM PLAMERE

Esthetic Multi Plasma Device
For Fibroblast Lifting Therapy

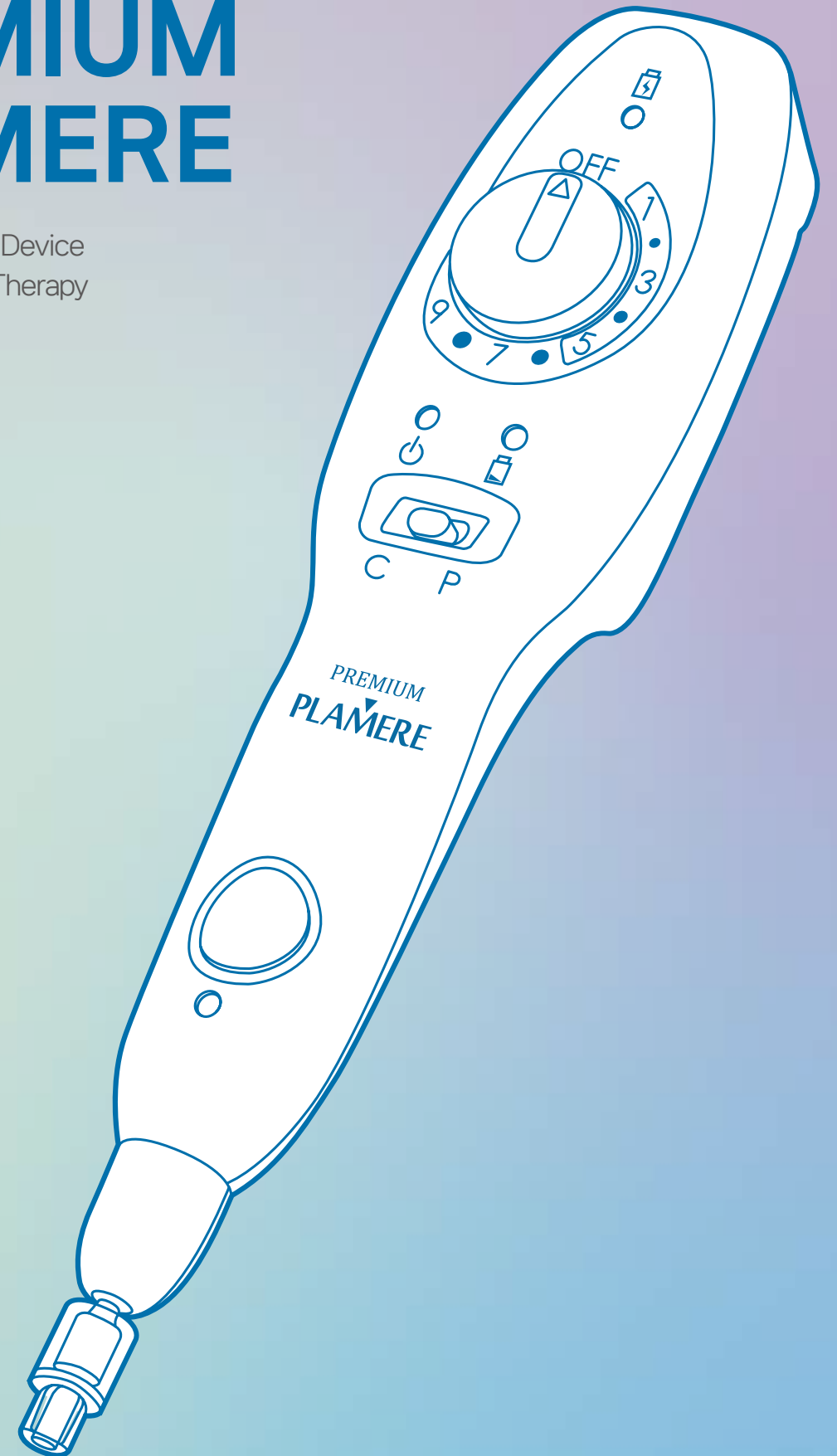

Product Name	PREMIUM PLAMERE
Model Name	GPM-2000
Output Frequency	40kHz
Output Level	1-9 (Analog Wheel Type)
Power	Internal Li-Polymer Rechargeable Battery, 3.7V, 500mAh
Power Consumption	Max, 2W
Charger	Micro USB
Mode	Continuous/Pulse
Usage Time	Max 3 Hours after full charging
Dimension & Weight	210 x 140 x 70mm (W x L x H) 450gm

GCS Co., Ltd.

Seoul office 8F, Magic Kingdom Bldg., 18, Teheran-ro 83-gil, Gangnam-gu, Seoul, 06169 Korea
Headquarter 3F, 33, Sagimakgol-ro 62beon-gil, Jungwon-gu, Seongnam-si, Gyeonggi-do, 13211, Korea
(Sangdaewon-dong, Center M)
Email sales@gcshcb.com

GCS
We make the healthy beauty together
Healthcare

PREMIUM PLAMERE

Traditional Plasma Pen for Fibroblast

Premium Plamere is the ultimate choice for discerning professional estheticians. With its innovative design featuring five distinct functional tips, this device unlocks a world of possibilities for an array of skin treatments. From plasma lifting to deep wrinkle removal, acne care, mesotherapy using fractional and permeating tips, Premium Plamere delivers exceptional versatility and astounding results.

Power Indicator

Indicates Power ON/OFF status.

Wheel Switch

Turn ON/OFF the power or adjust the plasma output intensity.

Battery Warning Light

Red is displayed when the remaining battery is low. If the warning light comes on, stop using it and charge it before using it.

Mode Selection Switch

Continuous (C) or Pulse (P) mode. If pulse mode is selected, the output will be in pulse form of 500Hz. When using the pulse mode, set the wheel switch 5 or less.

Shot Button

Operating button for plasma output.

Output Indicator



Indicates the state of plasma output.

Head Cap

Cap for fixing needle tip.

Collet (Inside head cap)

A part for inserting the needle tip as 3Φ standard size.

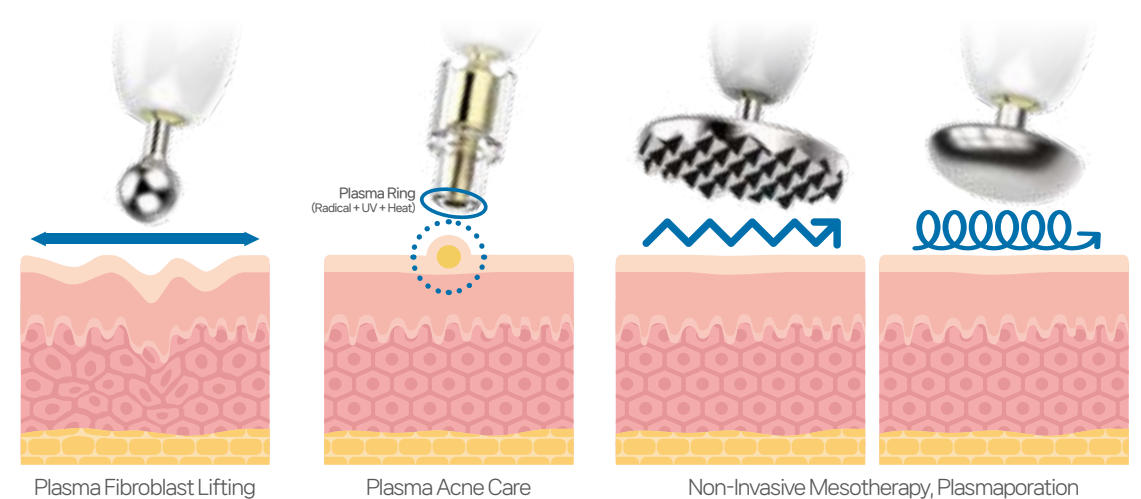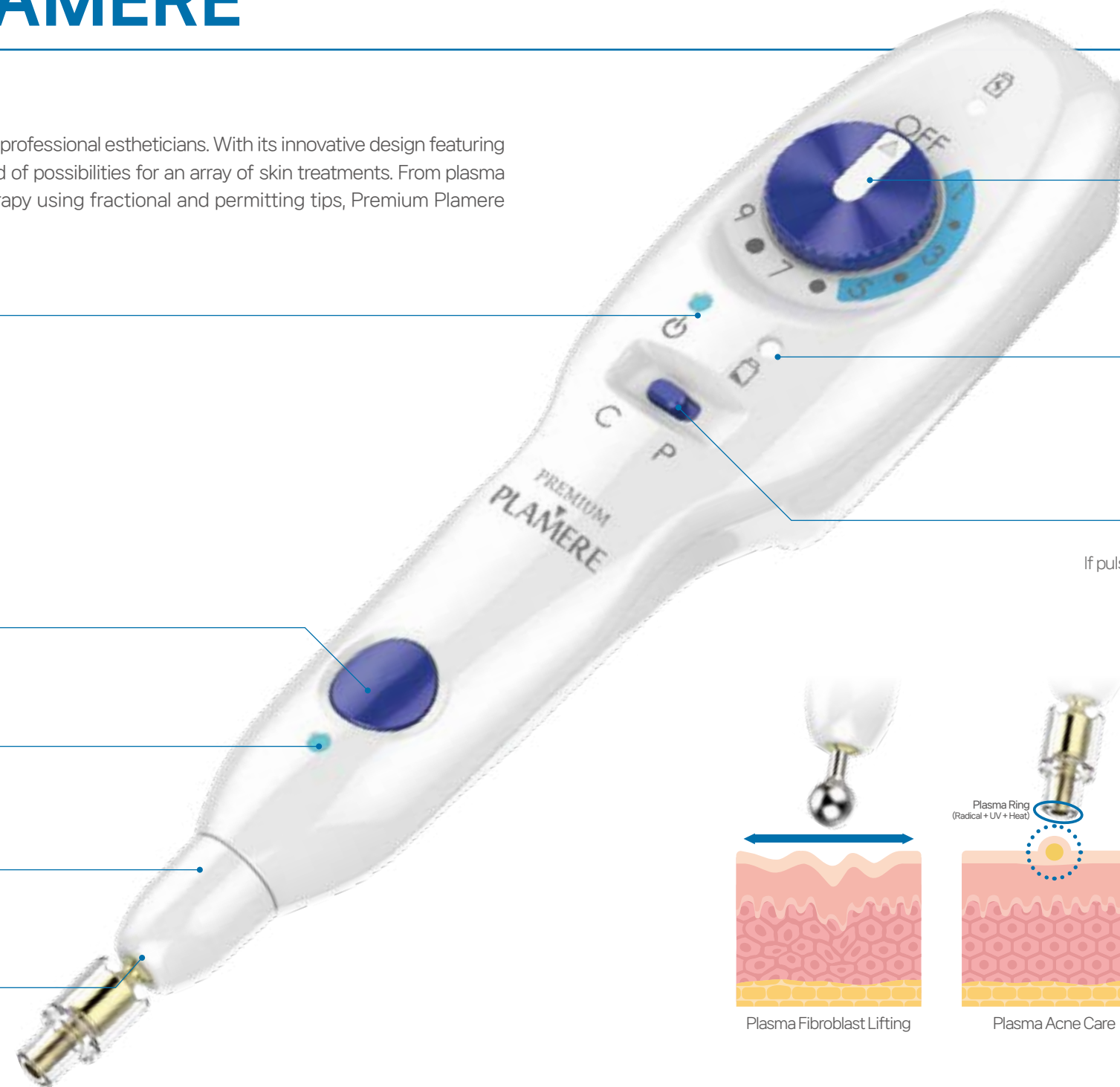

01 Plasma Fibroblast Lifting

Premium Plamere's Fibroblast lifting is a non-surgical aesthetic procedure that uses a high energy discharge to make a small wound in the skin. Our Plasma Arc stimulates increased activity of fibroblasts, collagen-producing cells in the connective tissue. So, it helps tighten and lifts the skin, face, and various parts of the body with instant results.

02 Plasma Acne Care : Maximum Acne Management

Premium Plamere is the only device that has a **plasma ring** for acne treatment, unlike other plasma pens. The plasma ring of Premium Plamere is highly effective in destroying Acne bacteria and reducing skin inflammation and preventing acne reoccurrence by balancing sebum production.

03 Non-Invasive Mesotherapy, Plasmaporation

Premium Plamere's tips help serum to be absorbed deep into the skin layers without invasive procedure by needle. Our fractional and permeating tips work safely and non-invasively transfer the substances for mesotherapy into your skin to treat a variety of problem areas.


04 World's First Detachable External Battery Pack

In a world-first innovation, **PREMIUM PLAMERE** comes with a detachable external battery pack. Enjoy extended usage by simply replacing the battery pack.

Ver 2.
Plasma Therapy Series

LEAF

Fusion Plasma Device
For Absorption Boosting & Lifting Therapy

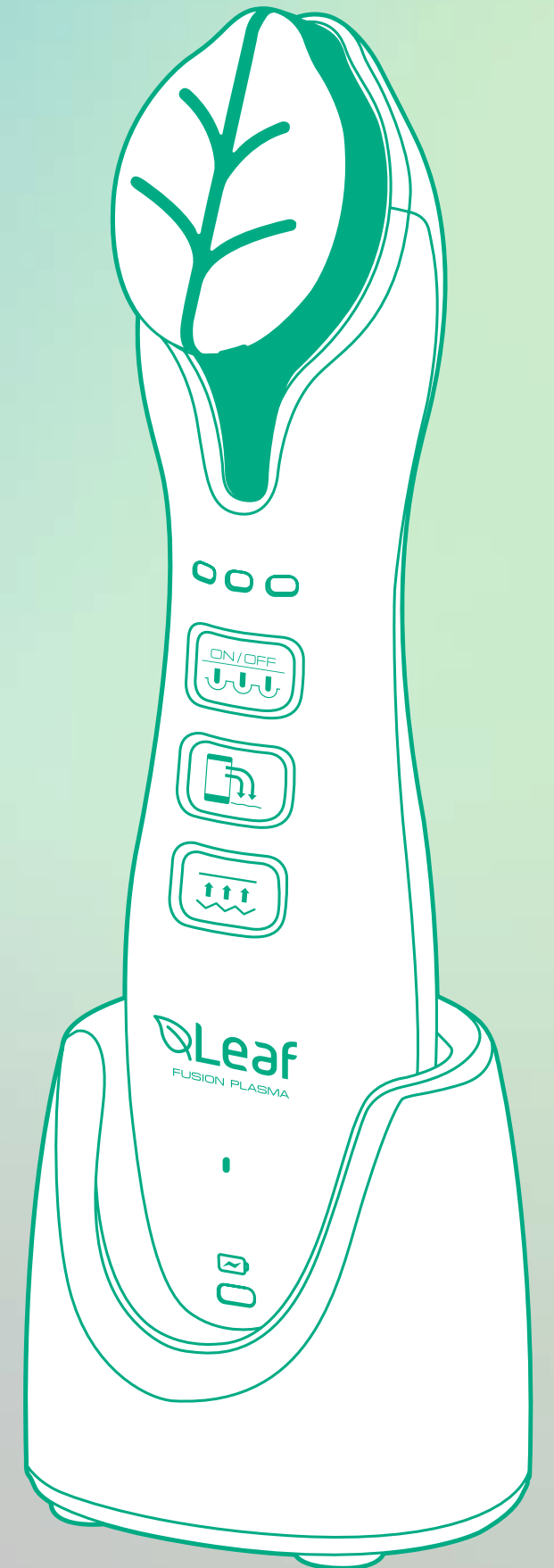

Product Name	LEAF FUSION PLASMA
Model Name	DFP-2000
Level Control	Level 1, Level 2, Level 3 [3 Levels]
Control Button	Fusion Plasma, Smart Boosting, Micro Current
Power	Built-in lithium polymer 3.7V 500mAh
Power Consumption	1W max
Usage Time	10 days when using 10 minutes per day [If fully charged]
Charging Type	Standard micro USB (Connect with computer USB port or mobile phone charger)
Product Composition Size	Main body, cradle, micro USB charging cable 42(W)x190(H)x40(D)(mm)
Weight(main body)	120g

GCS Co., Ltd.

Seoul office 8F, Magic Kingdom Bldg., 18, Teheran-ro 83-gil, Gangnam-gu, Seoul, 06169 Korea
Headquarter 3F, 33, Sagimakgol-ro 62beon-gil, Jungwon-gu, Seongnam-si, Gyeonggi-do, 13211, Korea
(Sangdaewon-dong, Center M)
Email sales@gcshcb.com


Ideal for Both Professional Beauticians and At-homecare Users

Revolutionize your self-care experience with Leaf, fusion plasma device designed to transform your skincare and well-being regimen. Packed with state-of-the-art features including plasma therapy, electroporation (EP), microcurrent stimulation, and IoT connectivity, Leaf is your all-in-one solution for radiant skin and wellness – all within the comfort of your own home.

Electroporation Excellence

Remarkable benefits of LEAF :

Effortlessly absorbs over 100 times more than applying it by hand.

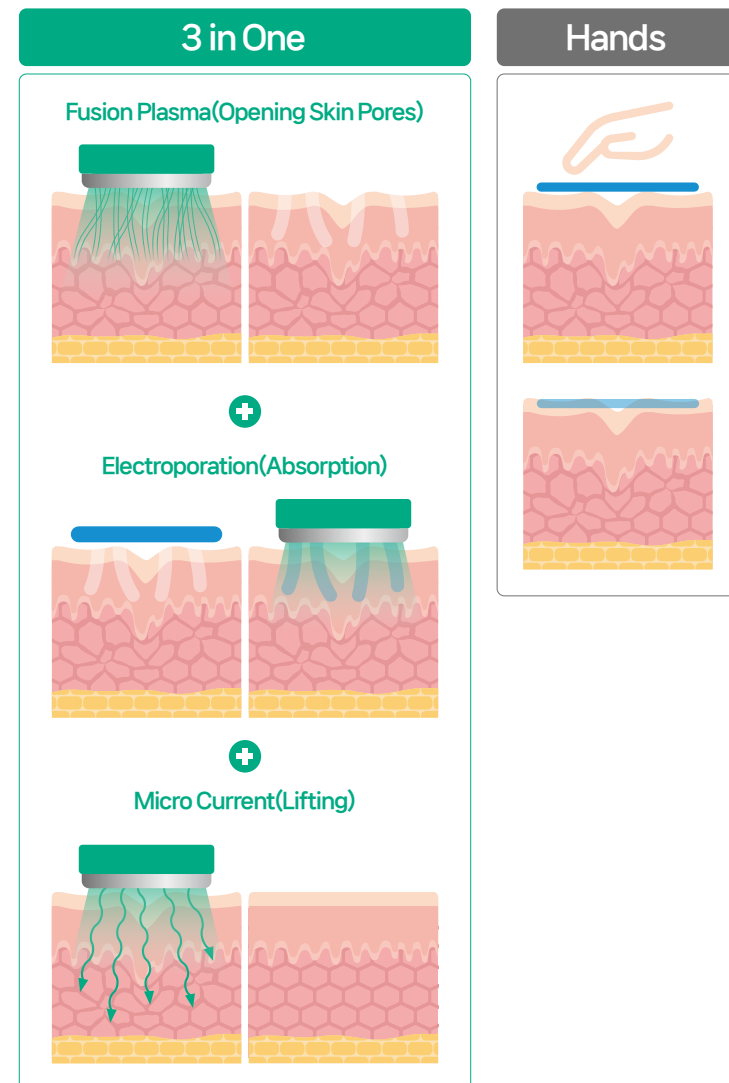

Head

Level Indicator

Main Body Power / Fusion Plasma Mode

Communication / Smart Boosting Mode

Voice Guide / Lifting Mode

Charging Cradle

Battery Status



IoT Connectivity :

Seamlessly integrate Leaf into your smart home ecosystem. Track your progress, receive personalized recommendations, and optimize your self-care routine through intuitive connectivity.

NO	Skin Type	Mode/Purpose	Head Color
1		Wrinkle	Blue
2		Soothing	Yellow
3	Sensitive Skin	Whitening	Purple
4		Moisturizing	Orange
5		Regeneration	Green
6		Wrinkle	Dark Blue
7		Soothing	Pink
8	Oily Skin	Whitening	Light Green
9		Moisturizing	Dark Blue
10		Regeneration	Light Green
11		Wrinkle	Dark Blue
12		Soothing	Purple
13	Dry Skin	Whitening	Orange
14		Moisturizing	Dark Blue
15		Regeneration	Purple
16		Wrinkle	Pink
17		Soothing	Light Green
18	Normal Skin	Whitening	Dark Blue
19		Moisturizing	Light Green
20		Regeneration	Dark Blue